

ENDOPULS 811

RADIAL SHOCKWAVE THERAPY



PASSION FOR MOTION



Partner for Life

Endopuls 811
the ideal unit for radial shockwave therapy,
suitable for every practice



ENDOPULS 811



The ideal unit for radial shockwave therapy, suitable for every practice

What is shockwave therapy?

Shockwave therapy is successfully used for the treatment of muscle/tendon tissue and tendon insertions.

Shockwaves are sound impulses with a lot of energy, which transfer as waves at depth. The energy released in the treatment area improves the metabolism and blood circulation, thereby stimulating the body's own repair mechanism.

Shockwave has many advantages. A few examples are:

- Short treatment time per session
- High percentage of successful therapies after relatively few treatments
- Renders medication and/or surgery superfluous

The effectiveness of shockwave therapy is sufficiently proven by scientific research.

What is the Endopuls 811?

The Endopuls 811 is a compact, and complete shockwave unit. This ideal shockwave therapy device is safe, comfortable and user-friendly, as well as cost-effective and space-saving, making it suitable for every practice.

Compact

The Endopuls 811 doesn't need a compressor, as it has an electromagnetic system. The electromagnetic system is small, light and mobile and requires little maintenance. The projectile in the ergonomic handpiece is electromagnetically activated.

Mobile

The Endopuls 811 is supplied in an attractive aluminium case with all the necessary accessories, such as various applicators and a foot switch. Easy to transport and everything at hand for a complete treatment.

Little maintenance, long lifetime

The Endopuls 811 needs very little maintenance. The ergonomic handpiece can be used for at least a 2 million waves (about 1000 treatments) before it needs revising.

The Endopuls 811, the complete unit

- The handpiece has an ergonomic design and is made of a heavy quality aluminium. Due to the handpiece's own weight less pressure is required by the therapist, it gives better intrinsic and vibration cushioning to protect the therapist, and a perfect dosage of energy transfer to the patient.
- Extremely easy to use: operated via the clearly displayed touchscreen.
- The energy is set precisely and displayed in mJoules, informing the therapist about the energy transferred to the tissue.



ENDOPULS 811

Endopuls 811: Radial shockwave therapy, mobile, long life, easy maintenance, clear touchscreen display



RADIAL AND ULNAIR EPICONDYLITIS



TENDINITIS OF THE SHOULDER /
SHOULDER PROBLEMS



CALCIFIC TENDINITIS OF THE
SHOULDER



PATELLAR TENDINITIS





PATELLAR TENDINOPATHY



ACHILLODYNIA



PLANTAR FASCIITIS



HEEL SPURS

ENDOPULS 811

Endopuls 811: Radial shockwave therapy, mobile, long life, easy maintenance, clear touchscreen display



MYOFASCIAL TRIGGER POINT
THERAPY E.G. NECK



MYOFASCIAL TRIGGER POINT
THERAPY E.G. BACK



DUPUYTRENS DISEASE



TROCHANTERIC BURSTITIS



PERIOSTITIS / SHIN SPLINT
(CONDITION AFTER OVERLOAD)



THUMB BASAL JOINT ARTHRITIS /
RHIZARTRITIS

ENDOPULS 811

Endopuls 811: Radial shockwave therapy, mobile, long life, easy maintenance, clear touchscreen display

ORDERING DATA

1650900 Endopuls 811

STANDARD ACCESSORIES

1650800 Treatment head
 1650808 Treatment head holder
 1650803 6 mm applicator point
 1650802 15 mm applicator point
 1650801 25 mm applicator point
 1650804 10 Silicon covers
 1650805 Massage lotion
 1650806 Foot switch
 1650807 Test template
 1650809 Transport case
 1650751 Operating manual
 xxxxxxx Power cable

OPTIONAL ACCESSORIES

1468950 EN-CAR U

TECHNICAL SPECIFICATIONS

Dosage / Energy settings: 60 mJ to 180 mJ in steps of 10 and 180 mJ - 185mJ
 (at 16 Hz maximum 120 mJ and at 22 Hz maximum 90 mJ).
 In comparison with a pressure-driven unit this is: 1 / 2 / 3 and 5 Bar.
 on SD cards

Save and update:
 Treatment head: Ergonomic, anodised aluminium housing.
 Dimensions device: length 35 cm, wide 25 cm, height 14 cm
 Dimensions handpiece: length 23 cm, 5 cm diameter
 Weight device: appr. 2,7 kg
 Weight handpiece: appr. 0.85 kg (with cable)
 Durability: minimum of 2,000,000 shots
 End-pieces: 6 x 25 mm and 15 x 25 mm diameter, with a durability of 150,000 shots minimum

Electrical connection: 100 – 240 VAC / 50-60Hz 2.5A
Approvals: IEC/ EN 60601-1



STANDARD ACCESSORIES ENDOPULS 811



1650800



1650808



1650803



1650802



1650801



1650804

OPTIONAL ACCESSORIES ENDOPULS 811



1650805



1650806




1650807



1650809



1650751

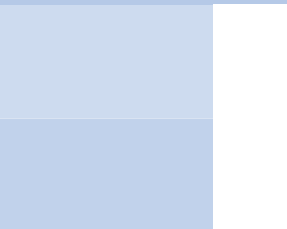

1468950

ENDOPULS 811



Endopuls 811: Radial shockwave therapy, mobile, long life, easy maintenance, clear touchscreen display

PASSION FOR MOTION

Enraf-Nonius B.V.
Vareseweg 127
3047 AT Rotterdam
P.O. Box 12080
3004 GB Rotterdam
The Netherlands
T +31 - (0)10 - 20 30 600
F +31 - (0)10 - 20 30 699
E info@enraf-nonius.nl
www.enraf-nonius.com

Technical modifications reserved

ENDOPULS 811